



Unison UK-C18 HT

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Sucralose
スクラロースの LC-MS/MS 分析

Standard 50 ng/mL

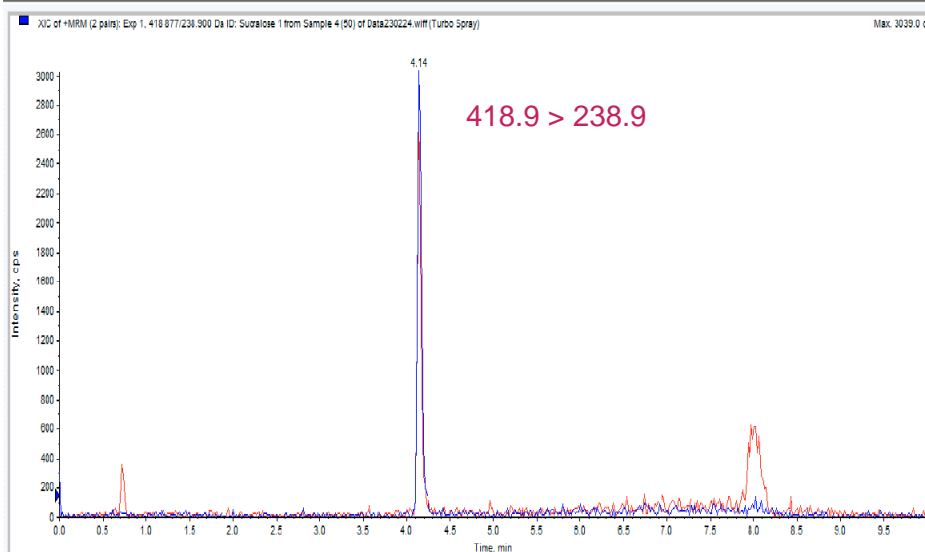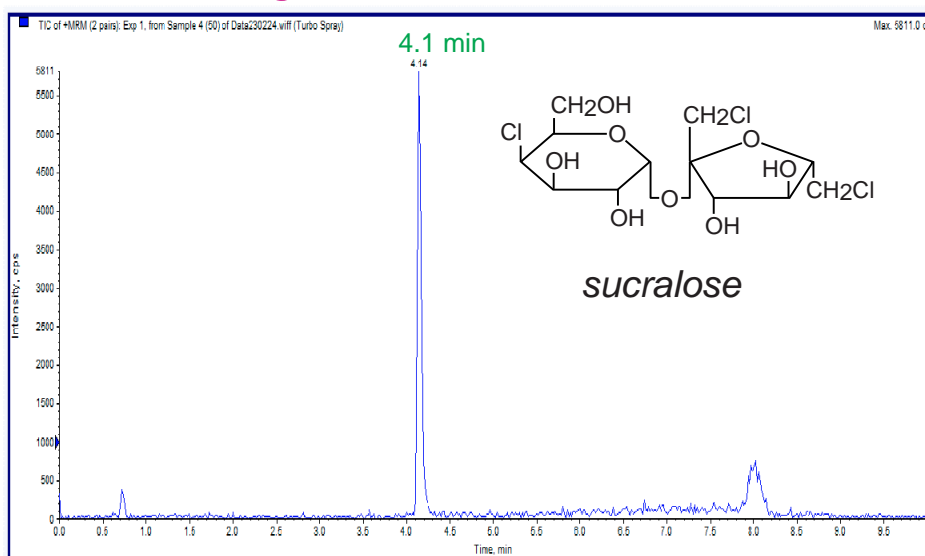

Unison UK-C18 HT, 100 x 2 mm
 A: 0.1% formic acid in water
 B: 0.1% formic acid in acetonitrile
 5%B (0-1.5min), 5-90%B (1.5-4min), 90%B (4-7min),
 90-5%B (7-7.1min), 5%B (7.1-10min)
 0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)
 2µL (50 ng/mL)

Courtesy of PEAKMAN SP, KOREA